


MILTEX® HI-LIGHT

ENDOMETRIAL SAMPLING SET IDEAL FOR SAMPLE COLLECTION OF CLEARLY DIFFERENTIATED ENDOMETRIAL TISSUE

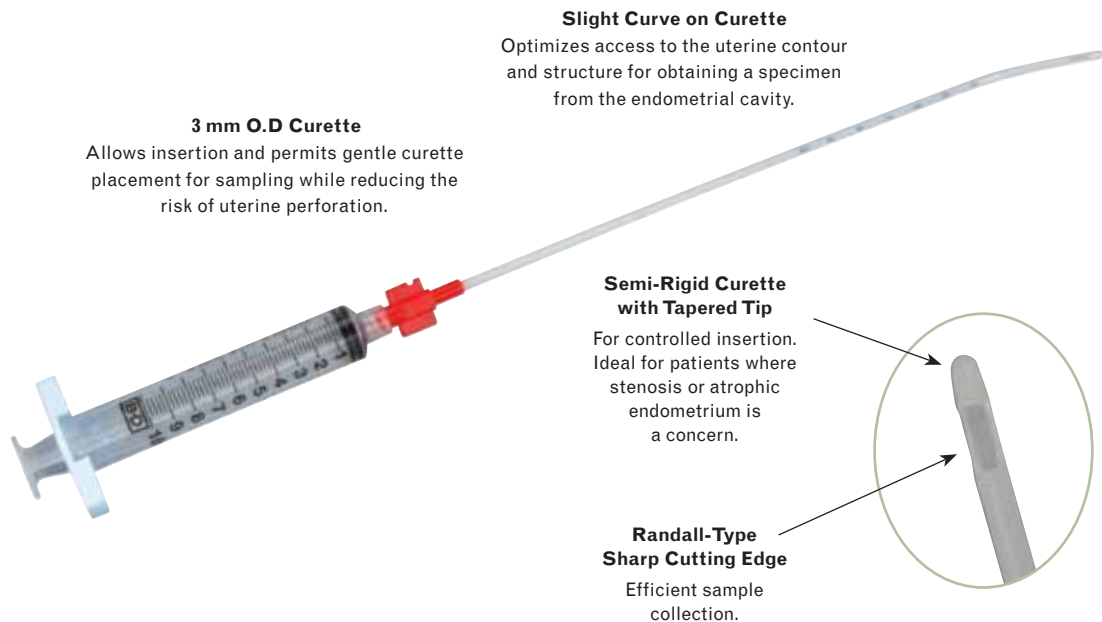

The Miltex® Endometrial Sampling Set is cost effective and accurate for screening of the endometrium. It is ideal for efficient sampling for both pre- and post-menopausal patients and affords reduced physical discomfort often associated with endometrial tissue collection.

FEATURES

- Set includes one each of Curette and Luer Lock Syringe
- 10 ml standard Twist-and-Lock Syringe
 - > Creates 21" Hg continuous negative pressure for optimized suction and sample capture
 - > Prevents backflow of sample into syringe
 - > Luer Lock fitting has winged flanges for ease of grasp and attachment
- Curette working length = 8-1/8" (20.64 cm)
- Disposable and latex-free
- Packaged 10 per box
- Sterile, 3-year shelf life

REF	Description
30-3015	Endometrial Sampling Set (includes Disposable Curette with a Twist-and-Lock Syringe), sterile, single-use, box of 10


See reverse side for additional Endometrial Sampling products.

Illustrations and content provide general description only and may be subject to change. For additional information, an in-office demonstration, or to place an order, please contact your local Miltex authorized Distributor or Miltex Customer Service at 1-866-854-8300, e-mail customerservice@miltex.com, or visit us at www.miltex.com.

Miltex and Miltex Redefining Excellence are registered trademarks of Miltex, Inc., a subsidiary of Integra LifeSciences Corporation. The Integra wave logo is a trademark of Integra Life Sciences Corporation. Copyright © 2010 Integra LifeSciences Corporation

ENDOMETRIAL BIOPSY CANNULA

TO OBTAIN SAMPLE FROM THE UTERINE MUCOSAL LINING

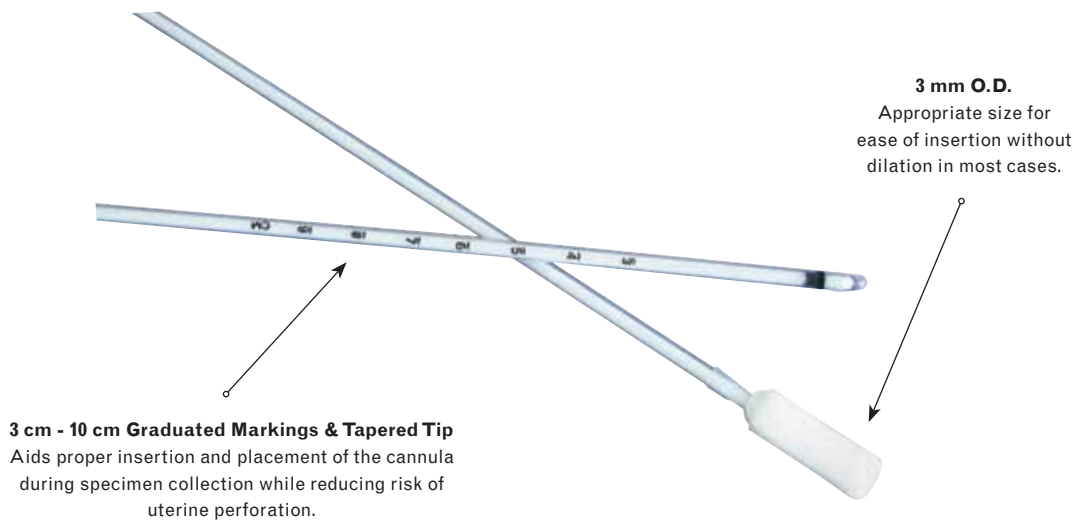

Endometrial sampling is the most widely used screening technique of the female patient population for a number of indications. The Miltex® SoftFlex™ and FirmFlex™ Endometrial Biopsy Cannula are ideal for accurate as well as economical in-office specimen collection of endometrial tissue.

TWO STYLES

SoftFlex™

- For standard/routine sampling
- When a more flexible cannula is desired

FirmFlex™

- For sampling when stenosis or atrophic endometrium is a concern
- For post menopausal women where stenosis of the cervix is a concern
- When a more rigid cannula is needed for controlled insertion

FEATURES

- One step sampling – no extra syringe or suction required
- Cannula and inner sleeve design generates sufficient vacuum for maximum sample yield without the need for external suction device
- 23.5 cm length allows for insertion depth to accommodate a range of patients
- Soft and Firm flexibility choices optimize formation to the natural uterine contour and structure for obtaining specimens from the endometrial cavity
- Disposable and latex-free
- Packaged 25 per box
- Sterile, 3-year shelf life

REF	Description
30-3010	SoftFlex™ Endometrial Biopsy Cannula, sterile, single-use, box of 25
30-3012	FirmFlex™ Endometrial Biopsy Cannula, sterile, single-use, box of 25