

Supplemental Phlebectomy Procedure Set

The following procedure set contains the most commonly utilized surgical instruments for Ambulatory Phlebectomy procedures. Integra™ Miltex® surgical instruments are manufactured from quality German stainless steel. These premium grade surgical instruments will provide desired clinical results during Phlebectomy procedures.



REF	Description	Qty
4-7	Knife Handle no. 3 fitting surgical blades Nos. 10 thru 15	1
7-10	Halsted Mosquito Forceps 5" (12.7 cm), curved, extra delicate	6
7-16	Halsted Mosquito Forceps 5" (12.7 cm), curved, 1 x 2 teeth	1
5-284	Metzenbaum Scissors 5" (12.7 cm), curved, delicate pattern	1
5-306	Iris Scissors 4 1/2" (11.4 cm), curved	1
18-SC-1474	Stevens Tenotomy Scissors 4 1/2" (11.4 cm), curved, sharp/sharp, SuperCut	1
8-7TC	Webster Needle Holder 4 3/4" (12.1 cm), smooth jaws, extra delicate, Tungsten Carbide	1
8-50TC	Crile-Wood Needle Holder 6" (15.2 cm), serrated jaws, Tungsten Carbide	1
6-121	Adson Tissue Forceps 4 3/4" (12.1 cm), 1 x 2 teeth and cross serrated platform	1
6-125	DeBakey-Adson Tissue Forceps 4 3/4" (12.1 cm), atraumatic tips 1.5 mm wide	1
11-76	Senn Retractor 6 3/8" (16.2 cm), blunt, double ended	1
3-200200	Miltex® Tray with Silicone Mats 7" x 11" (17.8 x 27.9 cm), double layer	1
3-800	Miltex® Instrument Care Kit System Spray Lube (3-700), Surgical Instrument Cleaner (3-720), Surgical Instrument Stain Remover (3-740), Instrument Care DVD (3-810), Nylon Brush and Stainless Steel Brush	1



integralife.com/integra-miltex

toll-free phone 866-854-8300

toll-free fax 866-854-8400

phone 717-840-9335

fax 717-840-9347

Please call Customer Service today for additional information, an in-office demonstration, to place an order or to receive other Integra Miltex product catalogs and information.

Illustrations and content provide general description only and may be subject to change

Integra and the Integra logo are trademarks of, and Miltex is a registered trademark of, Integra LifeSciences Corporation or its subsidiaries.

©2011 Integra LifeSciences Corporation. All rights reserved





IntegraTM Miltex[®]

Phlebectomy
Instruments

Crochet Phlebectomy Hooks

Crochet hooks are popular with clinicians due to the variety of sizes available. These hooks are a solid choice for exteriorizing the vein through the incision. Crochet patterns have a blunt, rounded tip and provide optimal elevation for all removal techniques. (German, stainless steel)

REF	Description	Size	Length
10403	Style 00	3.00 mm	6" (15.2 cm)
10404	Style 0	2.75 mm	6" (15.2 cm)
10405	Style 1	2.60 mm	6" (15.2 cm)
10406	Style 2	2.45 mm	6" (15.2 cm)
10407	Style 3	2.30 mm	6" (15.2 cm)
10408	Style 4	2.15 mm	6" (15.2 cm)
10409	Style 5	2.00 mm	6" (15.2 cm)
10410	Style 6	1.85 mm	6" (15.2 cm)
10411	Style 7	1.70 mm	6" (15.2 cm)
10412	Style 8	1.55 mm	6" (15.2 cm)
1040	Set, includes all 10 Styles #00 - 8	-	6" (15.2 cm)



10406



Mueller-Style Phlebectomy Hooks

Mueller-Style hooks have a blunt tip and flat handle which provide precise control during procedures. Mueller-style hooks are utilized to gently dissect the varicosities by undermining the vein laterally and perpendicularly. The vein is then exteriorized vertically using the point of the hook. Available in right and left handed patterns. (German, stainless steel)

REF	Description	Length
10321	Right #1	4-3/4" (12.1 cm)
10322	Right #2	5" (12.7 cm)
10323	Right #3	5-1/4" (13.3 cm)
10324	Right #4	5-1/4" (13.3 cm)
1032	Set, includes all 4 Right Hooks Sizes #1 - 4	-
10331	Left #1	4-3/4" (12.1 cm)
10332	Left #2	5" (12.7 cm)
10333	Left #3	5-1/4" (13.3 cm)
10334	Left #4	5-1/4" (13.3 cm)
1033	Set, includes all 4 Left Hooks Sizes #1 - 4	-




10321



Ramelet Phlebectomy Hooks

Ramelet hooks have a fine, sharp, angular hook and a short, round handle. The hook design facilitates dissection and provides a secure grip on the vein as well as moderate traction during vertical avulsion of the vein. This minimizes the damage to the surrounding tissues and lymphatics. The handle design allows for the instrument to be rolled between the fingers during dissection. This technique lessens the hand and wrist movements of the clinician during the procedure. Available in right and left handed patterns. (German, stainless steel)


REF	Description	Length
10351	Right #1	4" (10.2 cm)
10352	Right #2	4" (10.2 cm)
1035	Set, includes Right Hooks Sizes #1 - 2	-
10361	Left #1	4" (10.2 cm)
10362	Left #2	4" (10.2 cm)
1036	Set, includes Left Hooks Sizes #1 - 2	-

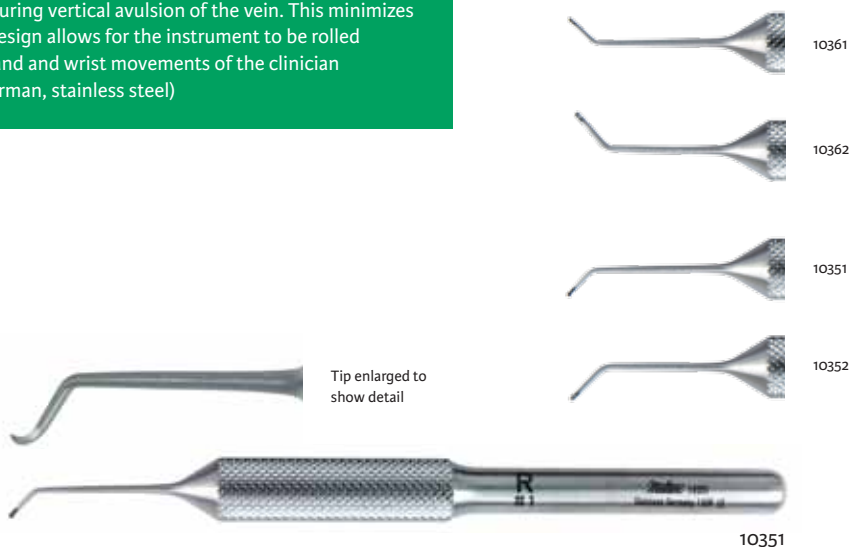

Oesch-Style Phlebectomy Hooks

Oesch-style hooks have a rectangular shaped handle for comfort and control during procedures. Small barbed tips are ideal for piercing the vein laterally for elevation. These fine tipped hooks permit usage through incisions as small as 1 millimeter to minimize damage to surrounding tissues and lymphatics. Three graduated tip sizes allow for the achievement of excellent cosmetic and clinical results. Available in right and left handed patterns. (German, stainless steel)

REF	Description	Length
10371	Right #1	6-1/2" (16.5 cm)
10372	Right #2	6-1/2" (16.5 cm)
10373	Right #3	6-1/2" (16.5 cm)
1037	Set, includes all 3 Right Hooks Sizes #1 - 3	-
10381	Left #1	6-1/2" (16.5 cm)
10382	Left #2	6-1/2" (16.5 cm)
10383	Left #3	6-1/2" (16.5 cm)
1038	Set, includes all 3 Left Hooks Sizes #1 - 3	-

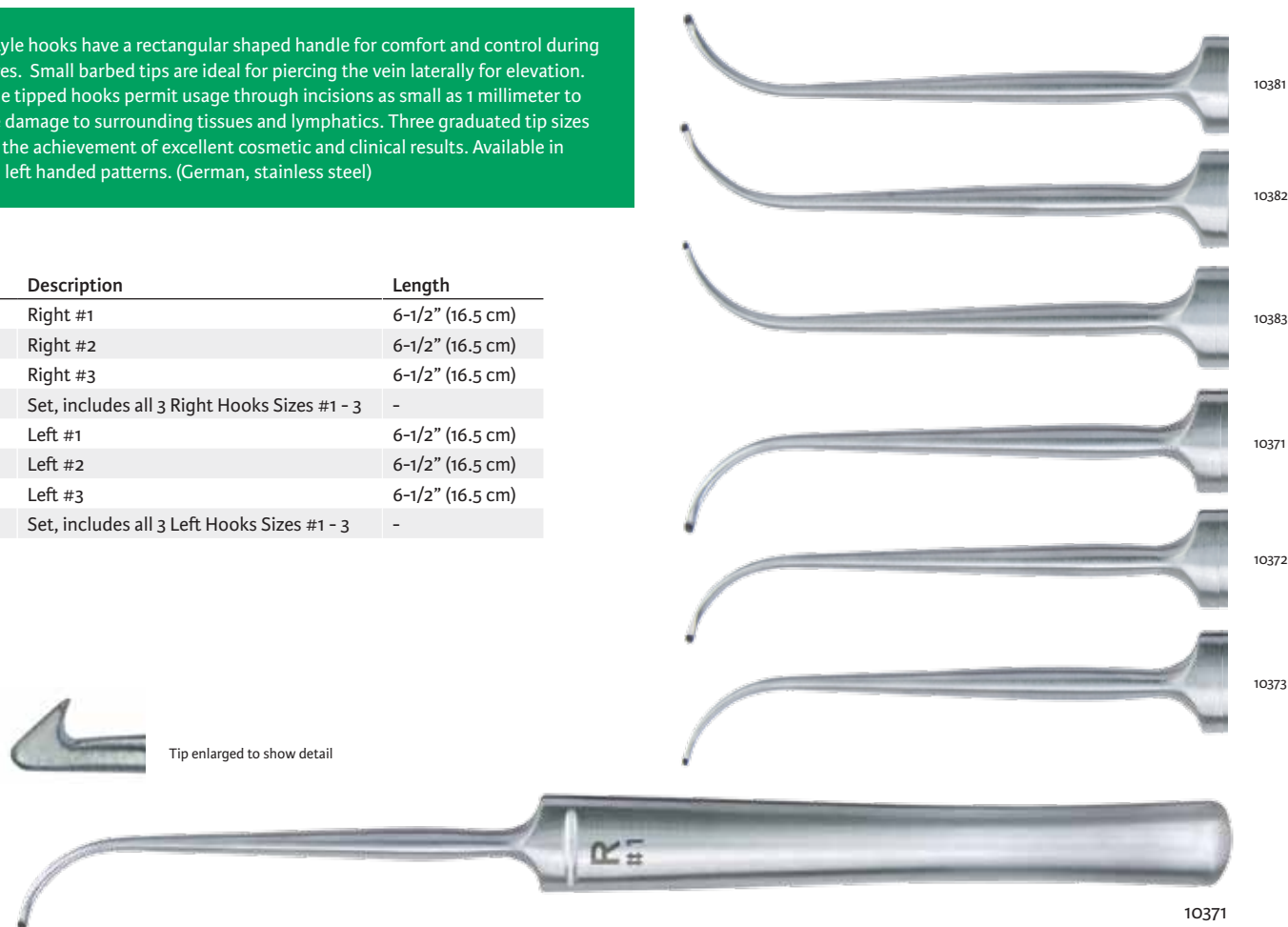

Varady Phlebectomy Extractors

Varady patterns are double ended with a variety of spatula and hook designs to meet all clinical needs. Spatula options help facilitate dissection and mobilization of the vein prior to elevation for extraction. Various hook designs and sizes are available to provide controlled elevation of the varix for exteriorization prior to proximal and distal avulsion. Specialization of Varady extractors provide optimal cosmetic and clinical results. (German, stainless steel)

REF	Description	Tips	Length
10305	Dissecting Spatula	1.9 mm and 2.9 mm width working ends	7" (17.8 cm)
10306	Large Hook w/Ball Tip and Spatula	2.8 mm wide hook and 1.9 mm wide spatula	7" (17.8 cm)
10307	Medium Crochet Style Hook and Spatula	1.8 mm wide hook and 1.7 mm wide spatula	7" (17.8 cm)
10308	Micro Spatula	0.7 mm and 1.2 mm width working ends	6-3/4" (17.1 cm)
10309	Micro Hook and Spatula	1.3 mm wide hook and 0.7 mm wide spatula	6-3/4" (17.1 cm)
10310	Small Hook and Spatula	1.5 mm wide hook and 1.1 mm wide spatula	6-3/4" (17.1 cm)
1030	Set, includes all 6 Extractors	-	-



10305



10306



10307



10308



10309



10310