

# MILTEX HI-LIGHT

## GRIP DEEP™ EXTRACTING FORCEPS



DEFDG1STD  
Extracting Forceps,  
Grip Deep 1 Standard

Grip Deep™ extracting forceps have been designed to allow subgingival access to the root while removing a tooth. The beaks are tapered, which allows them to penetrate deeper into the tooth socket and maximize the tooth-to-forceps contact. The beaks are also serrated which allows for better gripping of the tooth surface. This helps aid in atraumatic removal of the tooth and reduces the risk of fracturing the root structure, or cause trauma to the surrounding hard and soft tissues. The risk of crown fractures is also reduced.

The Miltex Grip Deep extracting forceps also have been designed with a take-apart feature that aids in the thorough cleaning and sterilization of the instruments. The forceps have a right and left side, which can be separated to assist in cleaning and sterilization of the hinge area. The Grip Deep extracting forceps are available in the 10 most popular patterns which include upper and lower molars, pre-molars, incisors and root tips.

Part Number	Description
DEFDG1STD	Extracting Forceps, Grip Deep 1 Standard
DEFDG99C	Extracting Forceps, Grip Deep 99C
DEFDG150	Extracting Forceps, Grip Deep 150
DEFDG150A	Extracting Forceps, Grip Deep 150A
DEFDG151	Extracting Forceps, Grip Deep 151
DEFDGMD3	Extracting Forceps, Grip Deep MD3
DEFDGMD4	Extracting Forceps, Grip Deep MD4
DEFDG210S	Extracting Forceps, Grip Deep 210S
DEFDG222	Extracting Forceps, Grip Deep 222
DEFDG65	Extracting Forceps, Grip Deep 65
DEFDG845	Extracting Forceps, Grip Deep 845

See reverse side for beak details

For additional information, an in-office demonstration, or to place an order, please contact your local Miltex, Inc. authorized Distributor or Miltex, Inc. Customer Service at 1-866-854-8300, email [customerservice@miltex.com](mailto:customerservice@miltex.com), or visit us at [www.miltex.com](http://www.miltex.com).



\*+H834DENTGDFCPHI0A\*

*Miltex*  
Redefining Excellence

phone 717 840-9335  
toll free 866 854-8300  
[www.miltex.com](http://www.miltex.com)

**DEFDG1STD**

Grip Deep 1 Standard  
*Upper Incisors & Cuspids*



**DEFDG99C**

Grip Deep 99C  
*Upper Incisors, Cuspids & Bicuspids*



**DEFDG150**

Grip Deep 150  
*Upper Incisors, Cuspids, Bicuspids & Roots; Universal*



**DEFDG150A**

Grip Deep 150A  
*Upper Incisors, Cuspids, Bicuspids; Universal; Parallel Beaks*




**DEFDGMD3**

Grip Deep MD3  
*Lower Incisors, Cuspids, Bicuspids and Roots, Mead; Universal*



**DEFDGMD4**

Grip Deep MD4  
*Lower 1st, 2nd and 3rd Molars, Mead; Universal*

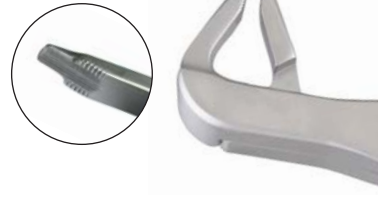

**DEFDG210S**

Grip Deep 210S  
*3rd Upper Molars; Universal*



**DEFDG222**

Grip Deep 222  
*3rd Lower Molars; Universal*



**DEFDG151**

Grip Deep 151  
*Lower Incisors, Cuspids, Bicuspids & Roots; Universal*



**DEFDG65**

Grip Deep 65  
*Upper Incisors, Roots & Fragments*



**DEFDG845**

Grip Deep 845  
*Lower Root Tips & Molars; Delicate*



**Serrated beaks aide in gripping and in the atraumatic removal of teeth.**



**Cup designed beaks to aide in root removal.**

